

Lab Table

Assembly Instructions



UVP, LLC
2066 W. 11th Street, Upland, CA 91786
Tel: (909) 946-3197 / (800) 452-6788
Fax: (909) 946-3597

Ultra-Violet Products Ltd.
Unit 1, Trinity Hall Farm Estate
Nuffield Road, Cambridge CB4 1TG UK
Tel: +44(0)1223-420022 Fax: +44(0)1223-420561

Web Site: www.uvp.com


81-0352-01 Rev STD

Assembly Instructions

It is essential that the screws are not tightened until all screws are inserted into their proper place. Once all screws are loosely in place, stand the Lab Table upright and tighten the screws.

1. Attach Vertical Legs A, F, G, and E to Middle Shelf D using the long, black 3/8 screws and nuts.
2. Attach Vertical Legs A, F, G, and E to Top Shelf B using the non-black 3/8 screws and nuts.
3. Attach Vertical Legs A, F, G, and E to Bottom Shelf C using the non-black 3/8 screws and nuts.

Schematic

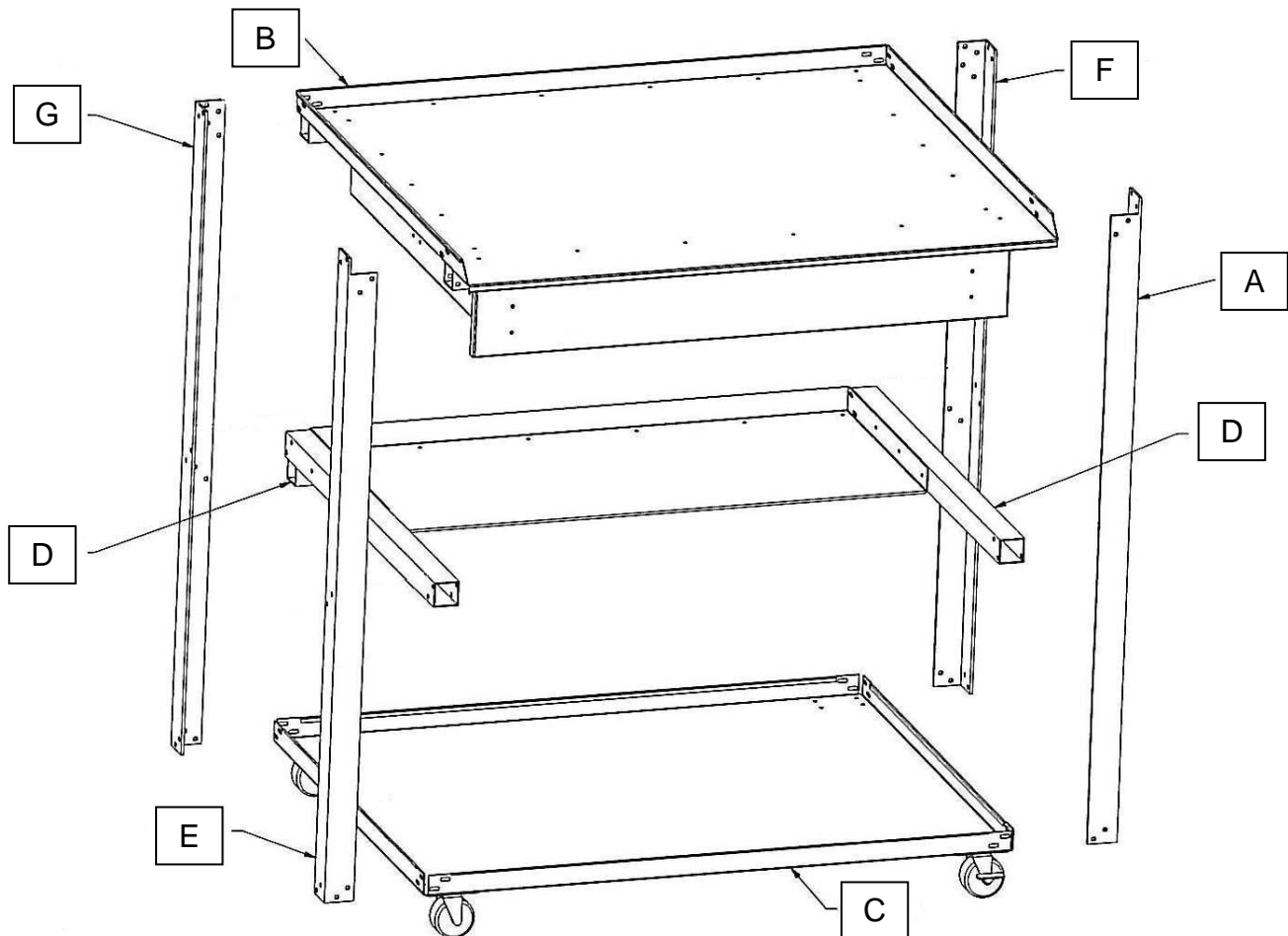

Corporate Headquarters: UVP, LLC
2066 W. 11th Street Upland, CA 91786 / Tel: (800) 452-6788
(909) 946-3197 / Fax: (909) 946-3597

European Operations: Ultra-Violet Products Ltd
Unit 1, Trinity Hall Farm Estate, Nuffield Road,
Cambridge CB4 1TG UK / Tel: +44(0)1223 420022
Fax: +44(0)1223 420156

For additional information on UVP's products, go to www.UVP.com